

THE “A” LIST

Anatomic Restoration:

The supinX™ Distal Biceps Repair System facilitates tendon reattachment to the native distal biceps footprint through either a single-head or double-head repair technique. By supporting anatomic tendon placement, the system is designed to recreate the distal biceps' natural biomechanics and function including forearm supination strength and rotational range of motion.

Anterior Single-Incision:


Employs the standard and widely accepted single incision, anterior approach supporting a consistent and reproducible technique for accessing the radial tuberosity and reattachment of the bicep tendon.

Atraumatic:

The supinX™ procedure utilizes a minimally invasive, 2mm bone-preserving tunnel to reattach the bicep tendon. The low-profile cortical button sits on the anterior aspect of the radial tuberosity, clear of neurovascular structures.


Automation:

The patented GOLDEN RETRIEVER™ guide with its self-capturing design provides precise targeting and streamlined suture management within a single workflow. Its automated retrieval of the suture from the non-visible, posterior aspect allows for restoration of the tendon to its native attachment, all from the anterior approach. The proprietary Button Threader automatically loads the suture button onto the suture strands, reducing fiddle, and simplifying the final step of securing the tendon construct to the bone. Together, these technologies are intended to support a controlled, efficient workflow and consistent suture management during distal biceps repair.

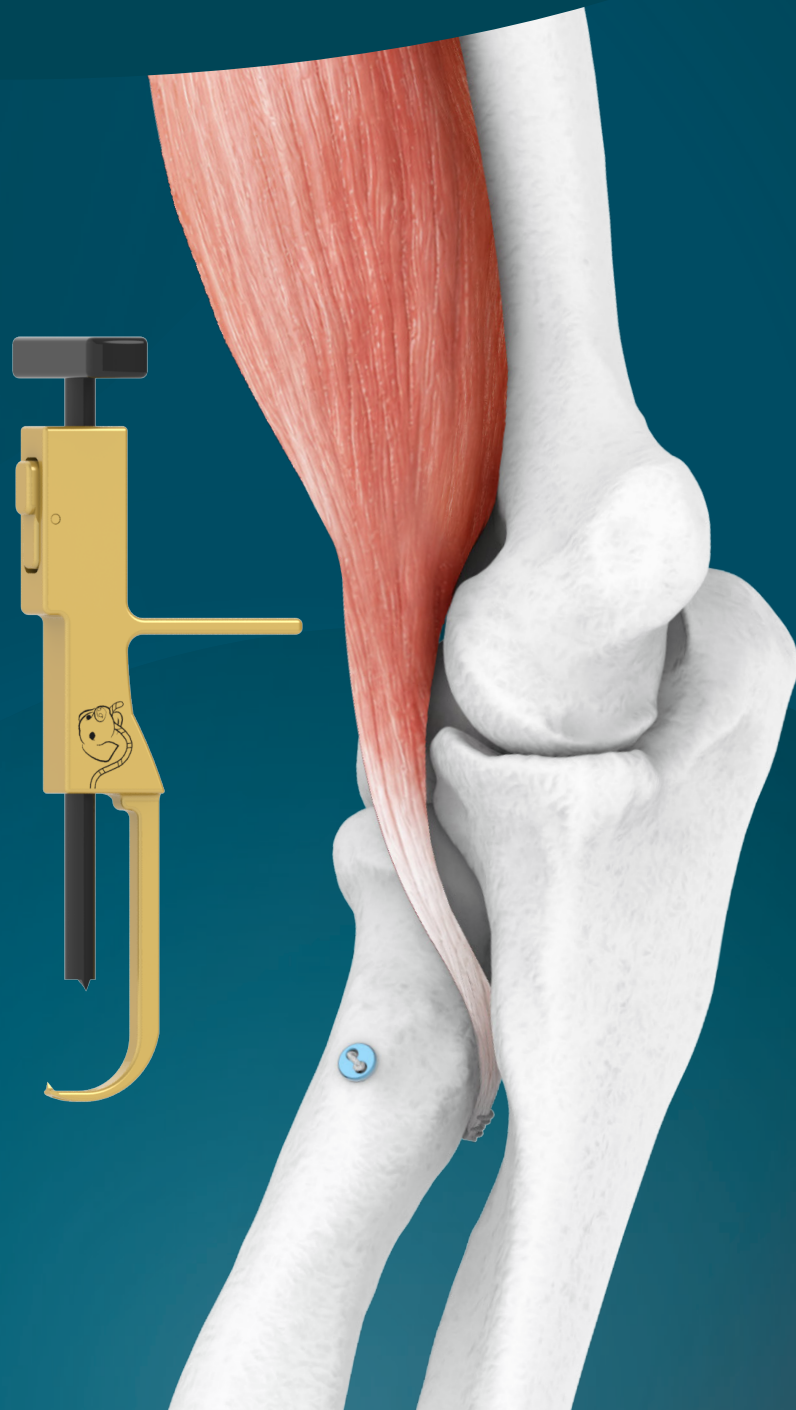

Enabling anatomic reattachment of the distal biceps tendon

supinX™ DISTAL BICEPS REPAIR SYSTEM


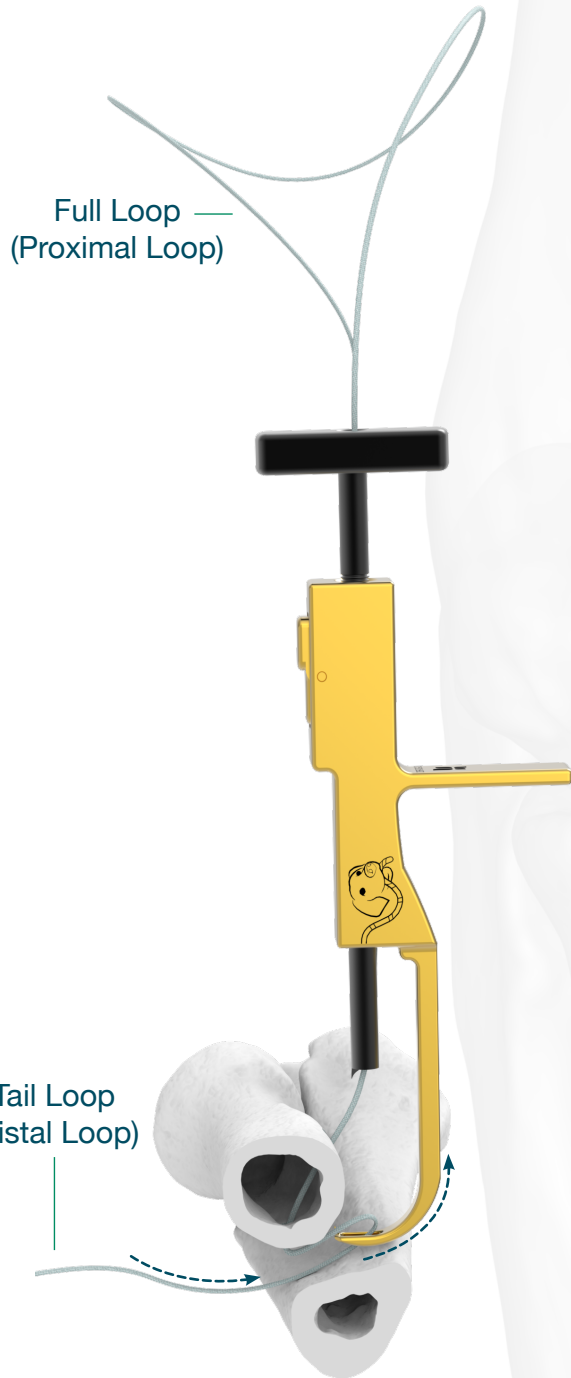


Anatomic re-attachment: Restores tendon to it's original insertion.

Anterior: Single incision approach.

Atraumatic: Minimal 2mm bone tunnel.

Automation: "Golden Retriever" guide retrieves

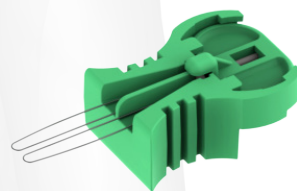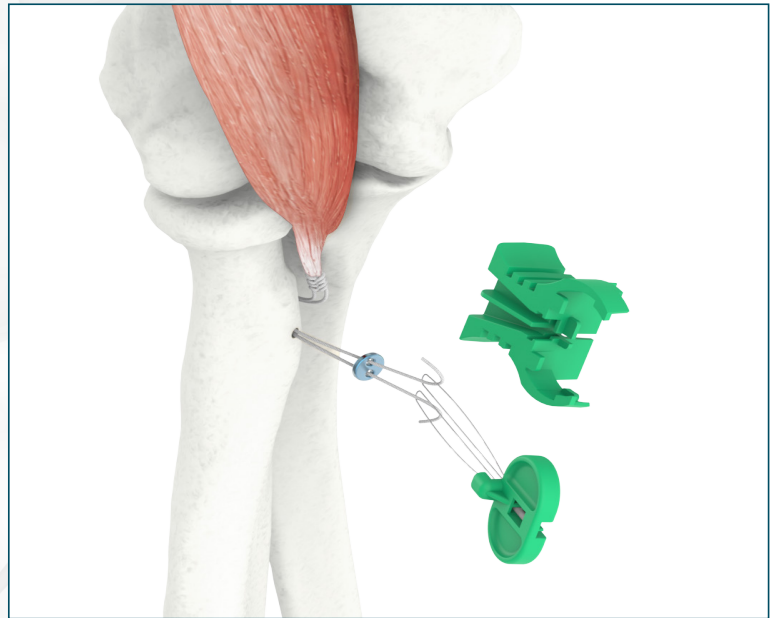

Button Threader:

Pre-loaded with Suture Button for convenient delivery to suture.

DISTAL BICEPS REPAIR SYSTEM CONSTRUCT

CATALOG #	DESCRIPTION
DBR-BLK-01 OR DBR-BLU-01	SupinX™ Distal Biceps Repair, Sterile Kit
DRLL-20124	Drill, 2.0mm x 124mm

THE KIT INCLUDES:

- #2 Suture (black or blue) – For tendon reattachment.
- Passing Needle – With Suture Loop.
- Green Button Threader – With preloaded Ti-Alloy Button.