



PROTEAN[®]

hand plating system




PERSONALIZED SOLUTIONS FOR INTERNAL FIXATION

The PROTEAN[®] Hand Plating System is designed for versatile stabilization and repair of metacarpal and phalangeal fractures, featuring specialized plates for intraarticular fractures of the metacarpal head and base, as well as multifragmentary fractures of the metacarpals and phalanges in children, adolescents and adults. The system features plates that are malleable in the vertical, horizontal and transverse planes that can be personalized on the back table or in-situ to fit the anatomic needs of each patient.



PROTEAN® hand plating system

PRE-CONTOURED STANDARD & FRACTURE SPECIFIC PLATES

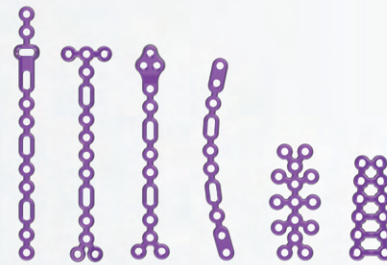

1.5 1.1 mm

2.0 1.4 mm

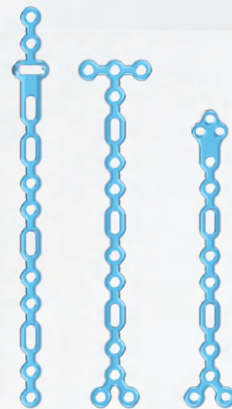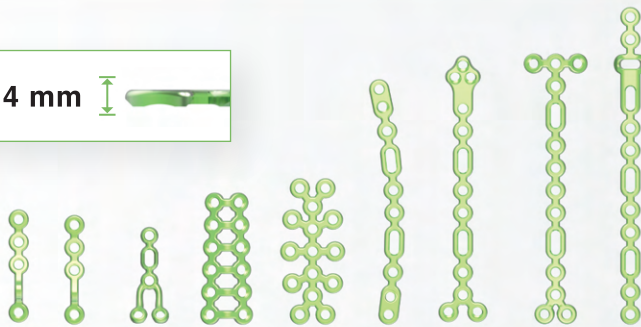

2.5 1.7 mm

NON-LOCKING & VARIABLE ANGLE LOCKING SCREWS

1.2



Non-Locking
CoCr

6-14mm

1.5



Non-Locking
Ti



Variable Angle
Locking CoCr

5-16, 18, 20 & 22mm

2.0



Non-Locking
Ti



Variable Angle
Locking CoCr

6-16, 18, 20 & 22mm

2.5



Non-Locking
Ti




Variable Angle
Locking CoCr

6-16, 18, 20 & 22mm



PROTEAN® hand plating system

PLATE PERSONALIZATION

The PROTEAN® Hand Plates provide exceptional versatility to address a broad range of anatomical and fracture reduction requirements. These plates can be precisely contoured and cut to the desired length using PROTEAN® Plate Bending Pliers, Bending Rods, or Cutting Pliers, providing patient and fracture specific personalization for each surgical case.


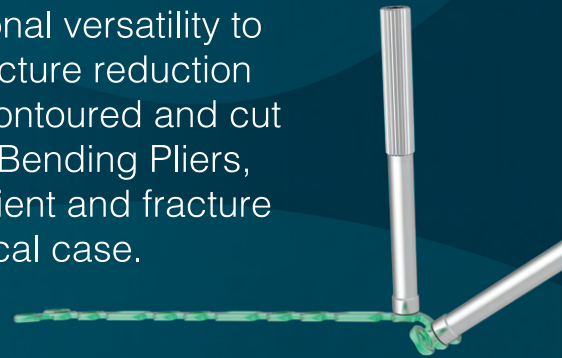


In plane bending
5° per node

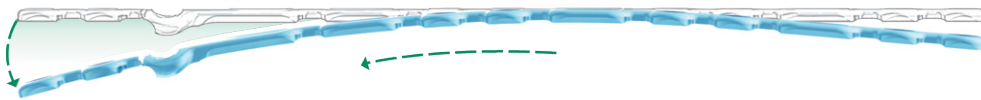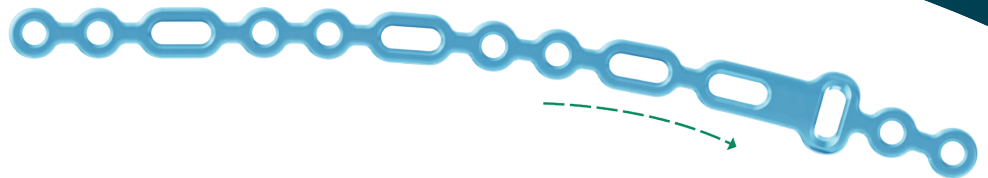

Out of plane bending
20° per node

Torquing
30° per node

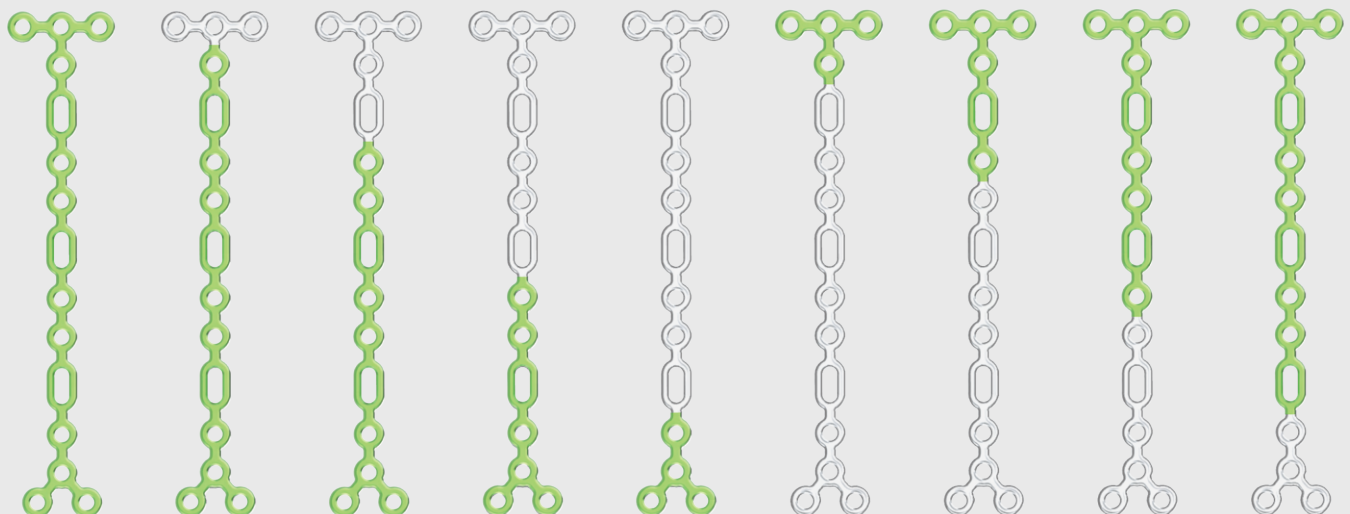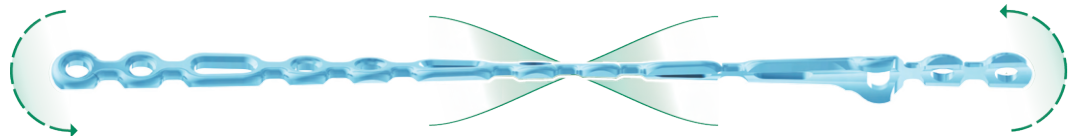

PROTEAN[®] hand plating system

UNIQUE PLATE + SCREW OPTIONS




Stability and Control from the Inside

Designed to preserve reduction and alignment while supporting early mobility.



IDEAL FOR DISTAL METACARPAL FRACTURES

Optimized for treating distal intra-articular metacarpal fractures, as well as extra-articular metacarpal fractures.

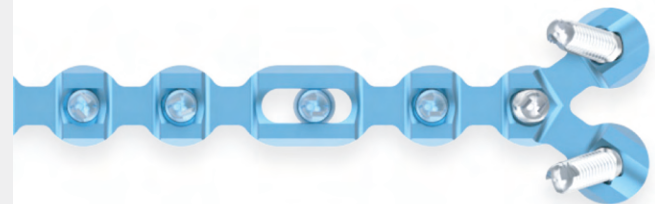


VERSATILE FOR PROXIMAL METACARPAL & PHALANGEAL FRACTURES

Designed to address proximal intra-articular metacarpal and phalangeal fractures, along with extra-articular fractures.



CURVED PLATE SOLUTION

Engineered for intra-articular fractures, as well as extra-articular fractures where a curved plate is required to place the plate from a lateral or medial approach.

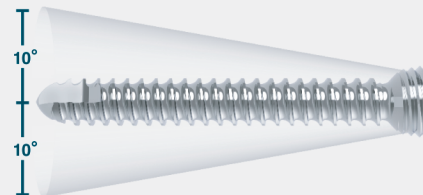

VARIABLE ANGLE LOCKING SCREWS FOR FLEXIBLE FIXATION

The advanced Variable Angle CoCr Locking Screws allow up to 10° off-axis angling in any direction.



PRECISION LAG SCREWS IN COCR

1.2mm CoCr lag screws are available for strength & durability in small bone applications.











PRECISION GUIDED CORRECTION AND FIXATION FOR CMC JOINT PRESERVATION

Created to preserve the CMC joint while alleviating pain and restoring thumb function.

PROTEAN® hand plating system

ONE TRAY, SIX MODULES

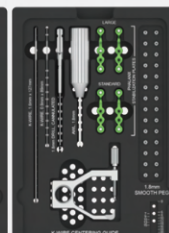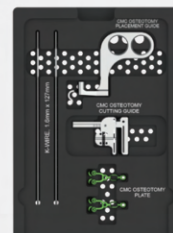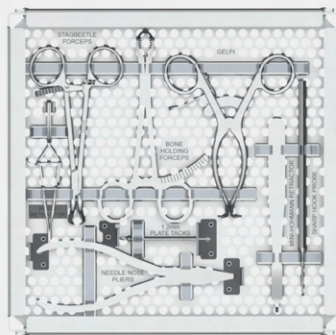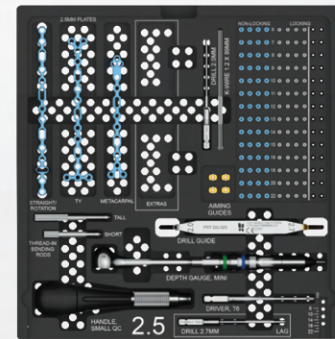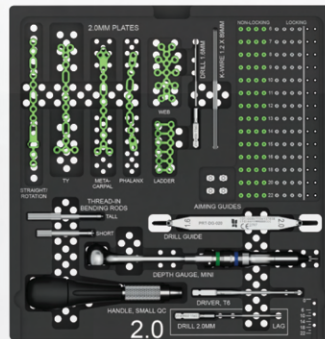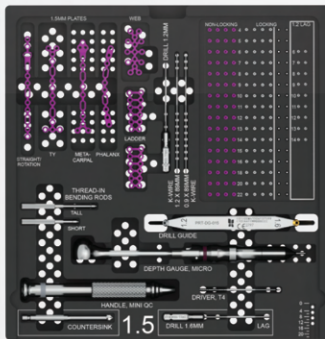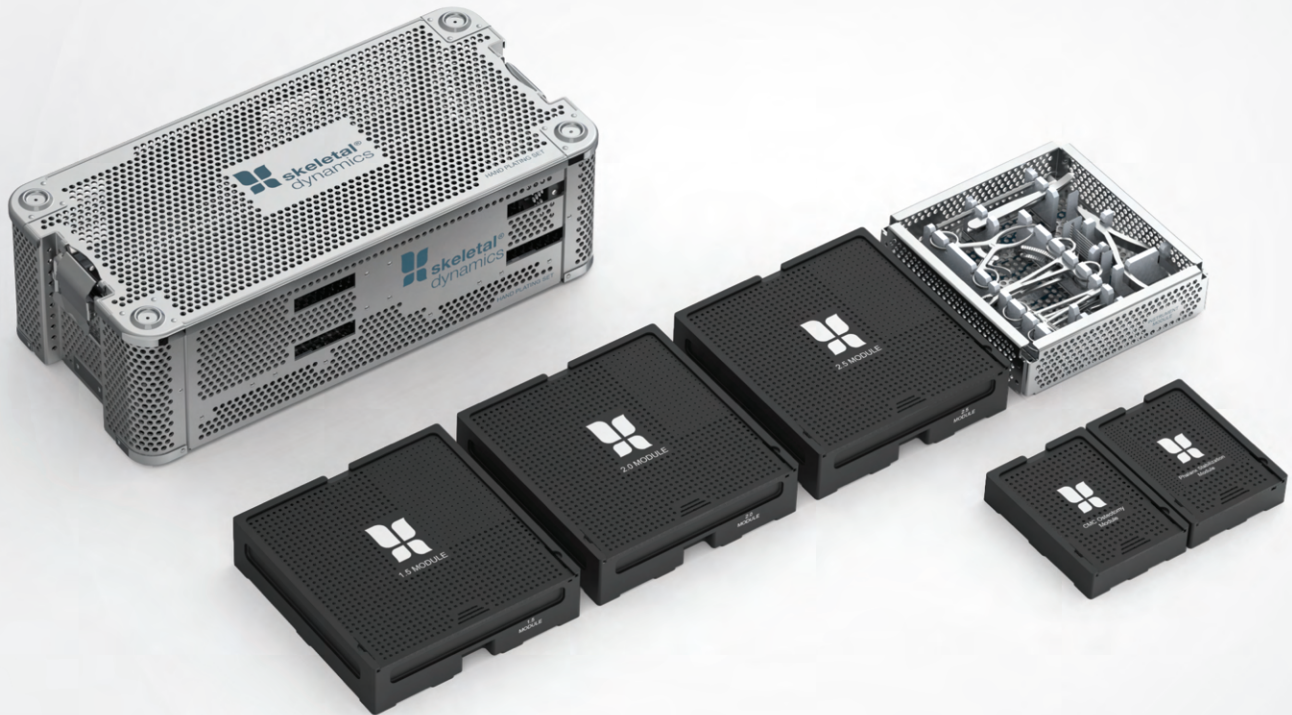

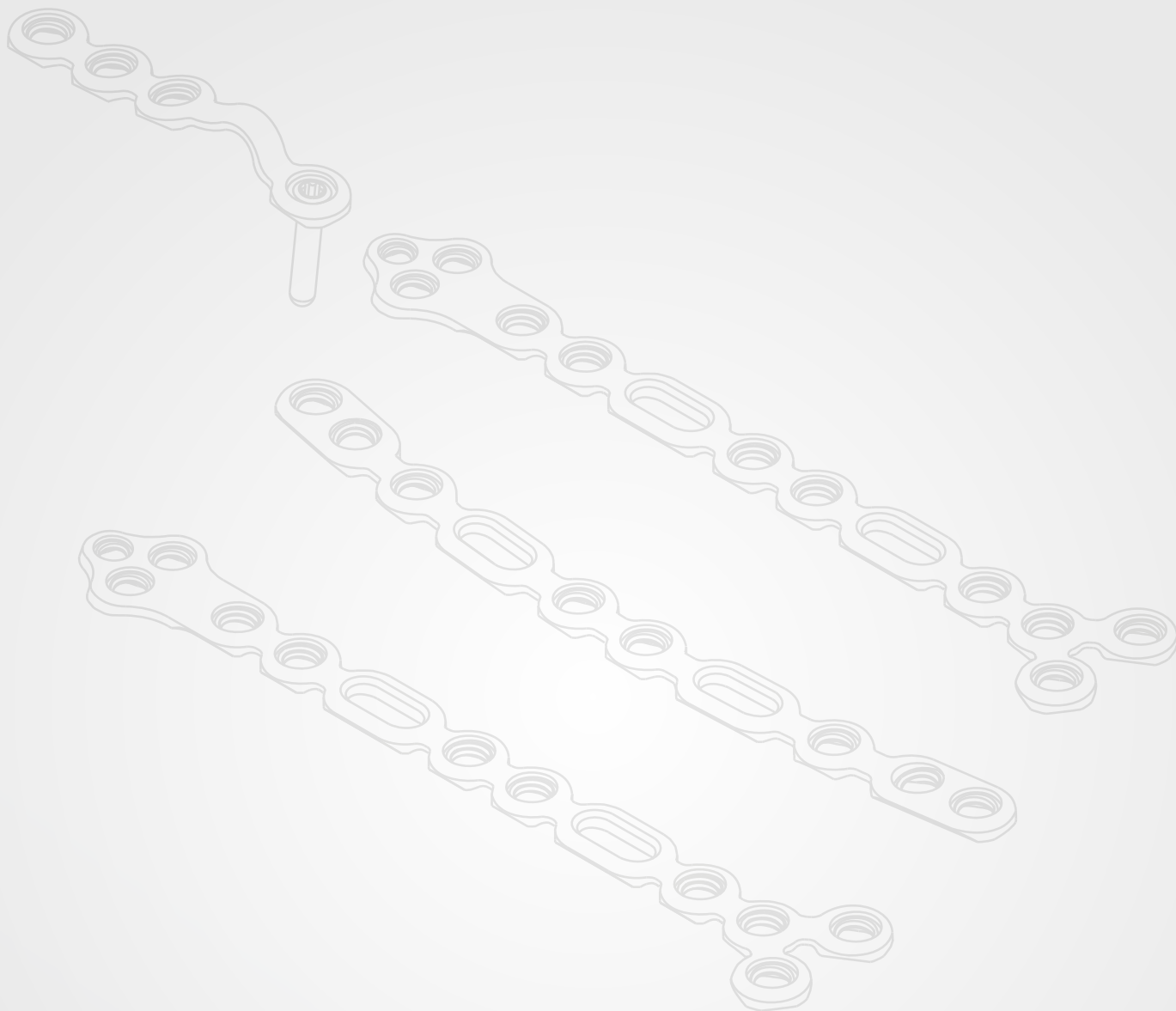

Scan the QR Code
Learn more about our
PROTEAN® Hand
Plating System

